

# ARCTIKO LR 100

true cooling specialists



## HIGH LEVEL OF SECURITY FOR YOUR SAMPLES, DUE TO THE GREAT COOLING PERFORMANCE WITH AIR FORCED COOLING.

The LR 100 will add a new dimension to your laboratory. The Freezer can be used as either a benchtop or undercounter model. The compact space-saving design and low noise level makes it a perfect storage solution for your laboratory.

### ► ATEX

This unit can be made in an ATEX version – approved according to EN 60079-15 category 3, zone 2. CE EX 11-3 G.

### ► FLEXIBLE STORAGE SOLUTIONS

Arctiko's LR range is a highly customizable allowing you to choose between a wide range of universal shelves and drawers.

This allows the user to organize their unit exactly as desired. Our universal shelf and drawer guides allow you to mix shelves and drawers for flexible storage.

### ► GREAT COOLING PERFORMANCE

Air forced cooling provides superior temperature uniformity and stability.

Forcing the airflow around the shelves and the stored items, ensures that all parts of the unit is equally cold and avoids any cold or hot spots.

Distributed by a powerful fan and open air flow areas along the sides, the air is always distributed evenly throughout the fully loaded unit.



LR 100\*



FORCED AIR COOLING



AVAILABLE IN ATEX



FILTERLESS CONSTRUCTION





FLEXIBLE SHELVES

\* The pictures shown are for reference only, actual product size may vary. Please read the dimensions of the product in the specifications.

The ARCTIKO LR 100 model supports following optional features:

- Modbus adapter for G214 controller
- Temperature chart recorder
- GSM alarm module
- Voltage converter
- Voltage stabilizer
- Stainless steel drawers
- Perforated shelf stainless steel
- Wire shelf
- LED Light

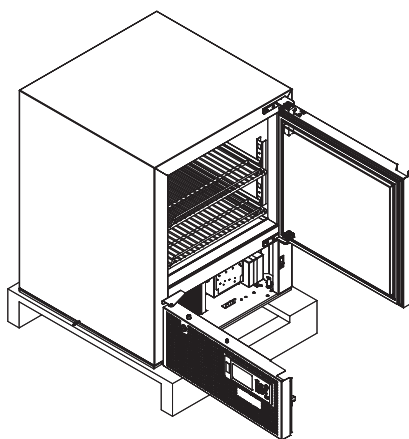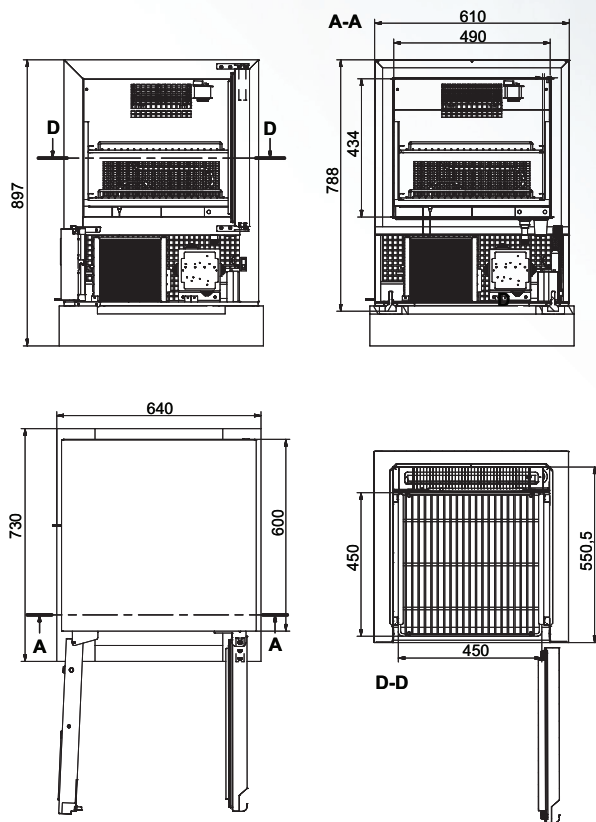

MODEL	LR 100	
Main specifications	Item code	DAI 0280
	Temperature range (°C)	+1 / +10
	Max. ambient temperature (°C)	43
	Capacity (L)	94
	Exterior dimensions WxDxH (mm)	610x659x790
	Interior dimensions WxDxH (mm)	490x440x434
	Weight (kg)	61
	Insulation (mm)	60
	Cooling technology	forced air
	Controller Model / probe type	G214 / PT1000
Power	Power supply (V)	230*
	Frequency (Hz)	50/60
	Power consumption (kWh/24h)	2.0
	AMP (A)	1.4
	Power (Watt)	206
	Horse power (HP)	1/8
Alarm functions	Alarm display as text - not codes	yes
	Visual / acoustic alarm	yes / yes
	Power failure alarm	yes
	Adjustable high / low temperature alarm	yes / yes
	Open door alarm	yes
	Probe failure alarm	yes
	Low battery alarm	yes
	Compressor failure alarm	yes
Battery backup for alarms (approx. hrs.)	72	
Data logging & external connection	Voltage free contact for remote alarm (e.g. GSM alarm module)	yes
	Electronic data logger	yes
	RS 485 port	yes
	USB port for software update & data download	yes
Temperature chart recorder	optional	
Display features	Display type	4,3" with touch buttons
	Number of optional reference probes	2
	Password protection for turning unit on / off	yes
	Temperature graph on display	yes
Other features	Automatic defrost	yes
	Antifreeze function	no
	Vacuum release port	no
	Access port for external probe etc. (int. diameter mm)	optional
	Light	no
	Lock with key	yes
	Castors (total pcs) / castors with brake (pcs)	2 rollers / 0
	Fixed feet (pcs) / levelling feet (pcs)	0 / 2
	Sub lids / inner doors (pcs)	no
	Reversible door(s)	no
	Fixed shelves (pcs) / adjustable shelves (pcs)	0 / 2
	Max. no. of shelves (pcs) / recommended max. no. of shelves (pcs)	5 / 3
Max. load per shelf (kg)	30	
Adjustable drawers (pcs)	optional	
Max. no. of adjustable drawers (pcs)	2	
Max. load per adjustable drawer (kg)	30	
Refrigeration & cabinet	Number of compressors	1
	Refrigerant	R290
	Interior cabinet material	stainless steel
	Exterior cabinet material	painting steel
	Colour (exterior cabinet)	white
	Foam type	polyurethane
VIP vacuum insulated panels	no	
Shipping	Shipping weight (kg)	83
	Shipping dimensions WxDxH (cm)	81x69x103
	Shipping volume (m³)	0,58

\*110V/60Hz available as option

